



Case number CAS-0022388810

Order number WO-31515749 | Visit Date 14/02/2025

Instrument Details

Instrument/Module cobas 6800

Serial Number 1376

Internal Instrument Name

Additional Details

PO

Customer Details

Customer Number 0052027574

Customer Name P'CHOMKLAO HOSPITAL
รพ.พระจอมเกล้า

Street Address 53 ถ.รถไฟ ต.คลองกระเซง

ZIP 76000 - Phetchaburi

Contact Name Supaporn Sae-aeaw

Contact Phone +66864684997

Service Activity Code Service Visit

Purpose of Visit

Instrument: C6800 SN 1376

Location: P'CHOMKLAO HOSPITAL รพ.พระจอมเกล้า

Issue: The instrument could not initialize, showing the message error with status of magazine C has changed to error. Suggested the customer to restart the instrument but the instrument cannot be back to standby and stuck on initializing status > Suggest to restart again, this time was back to standby but required PM to be performed. The customer performed the PM but it's stuck on initialization status.

Action: Tell me showed "The pusher drawer failed to initialize when it was not in power-off state. This could be due to the pusher-rubber-damper size at {3 CDI}"

The FSR suggested to close the reagent drawer loosely and the instrument was back to standby.

FSE was informed to further observe.

Performed Activities

clean & adjust motor step B/5/7

test drawer open / close ok

run test ok

Time Report



Category	Start	End	Hours	Billing Type
05913616001-Service labour time	14/02/2568 11:12	14/02/2568 16:00	4.80	Free of Charge

Signature or Stamp

Customer
Supaporn Sae-aeaw

Roche Support
Anurak Phaphai

Sign here

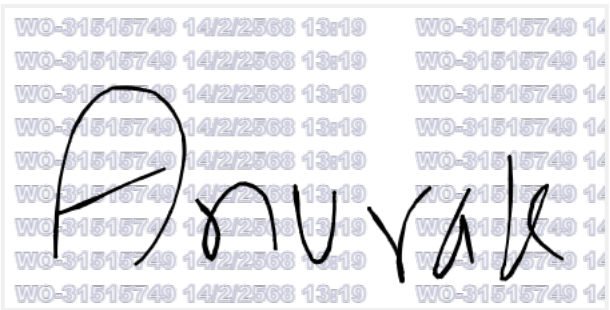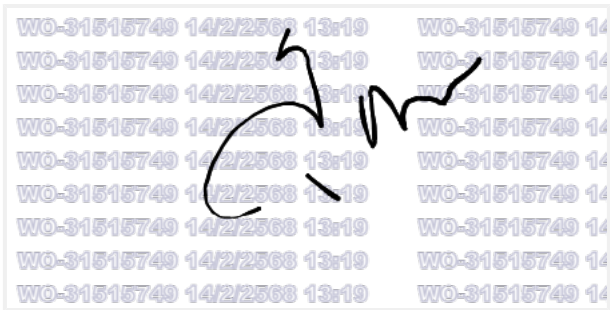

By signing this document

- Customer representative confirms they have reviewed the work carried out as described and accepts the system for use.
- Roche Support person confirms the above mentioned work and checks have been completed..
- .

Hotline : 1800-0193-99
Tel : 02 7912200
Line Id : @rcscthailand

Roche Customer Support Center