



Case number CAS-0023464304

Order number WO-33049743 | Visit Date 23/06/2025

Instrument Details

Instrument/Module cobas 6800

Serial Number 1376

Internal Instrument Name

Additional Details

PO

Customer Details

Customer Number 0052027574

Customer Name P'CHOMKLAO HOSPITAL
รพ.พระจอมเกล้า

Street Address 53 ถ.รถไฟ ต.คลองกระเซง

ZIP 76000 - Phetchaburi

Contact Name Supaporn Sae-aeaw

Contact Phone +66864684997

Service Activity Code Remote Support

Purpose of Visit

Related instrument : cobas 6800 s/n 1376 @P'CHOMKLAO HOSPITAL รพ.พระจอมเกล้า

Issue : A customer was complaint's about Transfer module was error after periodic maintenance. Error message was about reagent cassette drawer error. No error message code 1.3, PR was not detected any hardware issue.

ASA Tellme : Auto upload PR was detected :

Reagent Storage Sensor Anomaly 108

Anomalies 107/108/109 on reagent storage are triggered when a light bridge sensor pulls down the voltage below 24V.

ErrorNo 174 will subsequently occur on the next move.

Action : Restart instrument it was back to stand by. Inform FSE to RV to more inspection.

Event on 23/06/2025

Performed Activities

Related instrument : cobas 6800 s/n 1376 @P'CHOMKLAO HOSPITAL รพ.พระจอมเกล้า

Issue : A customer was complaint's about Transfer module was error after periodic maintenance. Error message was about reagent cassette drawer error. No error message code 1.3, PR was not detected any hardware issue.

ASA Tellme : Auto upload PR was detected :

Reagent Storage Sensor Anomaly 108

Anomalies 107/108/109 on reagent storage are triggered when a light bridge sensor pulls down the voltage below 24V.

ErrorNo 174 will subsequently occur on the next move.

Action : Restart instrument it was back to stand by. Inform FSE to RV to more inspection.

Event on 23/06/2025

Time Report


Category	Start	End	Hours	Billing Type
05913616001-Service labour time	23/06/2568 11:00	23/06/2568 12:00	1	Free of Charge

Signature or Stamp

Customer
Supaporn Sae-aeaw

Roche Support
Kasama Rakmark

Sign here

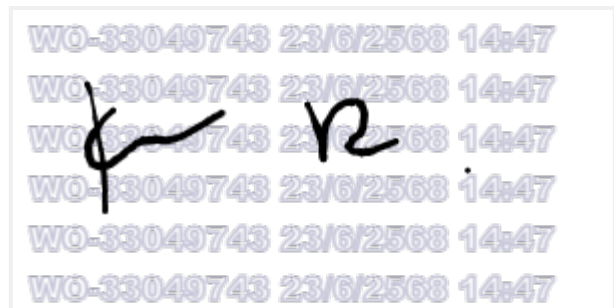

By signing this document

- Customer representative confirms they have reviewed the work carried out as described and accepts the system for use.
- Roche Support person confirms the above mentioned work and checks have been completed..
- .

Hotline : 1800-0193-99
Tel : 02 7912200
Line Id : @rcscthailand

Roche Customer Support Center