



Case number CAS-0023475192

Order number WO-33065337 | Visit Date 24/06/2025

Instrument Details

Instrument/Module cobas 6800

Serial Number 1376

Internal Instrument Name

Additional Details

PO

Customer Details

Customer Number 0052027574

Customer Name P'CHOMKLAO HOSPITAL
รพ.พระจอมเกล้า

Street Address 53 ถ.รถไฟ ต.คลองกระเซง

ZIP 76000 - Phetchaburi

Contact Name Supaporn Sae-aeaw

Contact Phone +66864684997

Service Activity Code Remote Support

Purpose of Visit

Instrument: cobas 6800 S/N 1376 P'CHOMKLAO HOSPITAL รพ.พระจอมเกล้า

Issue: The customer alleged transfer module has changed to error during performed cobas HPV run.

Action:

- Cleanup the analyzer >> standby
- check PR by using Tellme

Rootcauses

Detected Underaspiration:

The instrument detected insufficient aspirated sample volume . The test result was invalidated with flag U06T.

Sample transfer result - A "Failed" reusult indicates which head detected the underaspiration

Head A: Succeeded

Head B: Failed

Head C: Unused

Head D: Unused

Affected run batch id: 451ca403-b450-f011-80f8-000060e897f4

This is the first time that the customer found this issue (U06T Head B: Failed)

The customer will repeat run on 25 Jun 25. Then transfer case to FAS (Ms. Supatsa) for further investigation.

Performed Activities

Rootcauses

Detected Underaspiration:

The instrument detected insufficient aspirated sample volume . The test result was invalidated with flag U06T.

Sample transfer result - A "Failed" reusult indicates which head detected the underaspiration

Head A: Succeeded

Head B: Failed

Head C: Unused

Head D: Unused

Affected run batch id: 451ca403-b450-f011-80f8-000060e897f4

This is the first time that the customer found this issue (U06T Head B: Failed)

Time Report


Category	Start	End	Hours	Billing Type
05913616001-Service labour time	24/06/2568 14:30	24/06/2568 15:30	1	Free of Charge

Signature or Stamp

Customer
Supaporn Sae-aeaw

Roche Support
Monchanok Boonyahong

Sign here

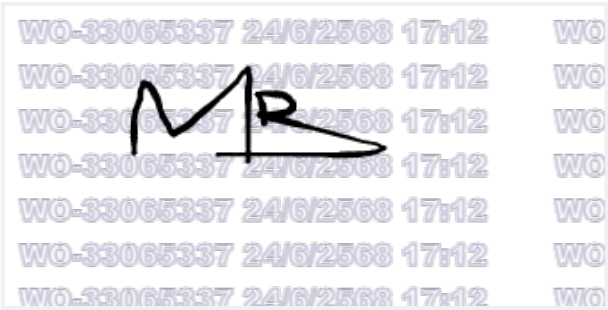

By signing this document
- Customer representative confirms they have reviewed the work carried out as described and accepts the system for use.
- Roche Support person confirms the above mentioned work and checks have been completed..
- .

Hotline : 1800-0193-99
Tel : 02 7912200
Line Id : @rcscthailand

Roche Customer Support Center