



Case number CAS-0023987933

Order number WO-33871065 | Visit Date 08/26/2025

Instrument Details

Instrument/Module cobas 6800

Serial Number 1376

Internal Instrument Name

Additional Details

PO

Customer Details

Customer Number 0052027574

Customer Name P'CHOMKLAO HOSPITAL
รพ.พระจอมเกล้า

Street Address 53 ถ.รถไฟ ต.คลองกระเซง

ZIP 76000 - Phetchaburi

Contact Name Supaporn Sae-aeaw

Contact Phone +66864684997

Service Activity Code Service Visit

Purpose of Visit

Instrument: S/N 1376 - cobas 6800

Location: P'CHOMKLAO HOSPITAL รพ.พระจอมเกล้า

Issue:

The processing module error occurred during an HPV run

Tellme:

- "Transfer Carrier failed"

During the transfer of a consumable or reagent cassette, an error occurred. The involved light barriers did not detect the correct states. A problem occurred in the ProcessCellA during the transfer of a consumable from 1 1 4002 13 7014 (TransferPosition) to 11 1 4002 13 7014 (Processing park1) with the reason TargetPositionOccupied.

- TargetPositionOccupied: Before the transfer started, the destination light barrier was already triggered when the consumable/reagent was still on the source position.

Action:

- Manual clean up, shut down the instrument to remove consumables
- Manual clean up again but did not standby > restart > passed
- Performed instrument teaching
- Informed K.Sompop for further investigation.

Affected test: 55 samples + 2 controls

Performed Activities

- Adjust Light Barrier Liquid Waste station 1
- Test Status Light Barrier
- initial pass
- standby run HPV test



Time Report

Category	Start	End	Hours	Billing Type
05913616001-Service labour time	08/26/2025 10:01 AM	08/26/2025 1:00 PM	2.98	Free of Charge

Signature or Stamp

Customer Supaporn Sae-aeaw	Roche Support Sompop Bangyeekhan
-------------------------------	-------------------------------------

Sign here

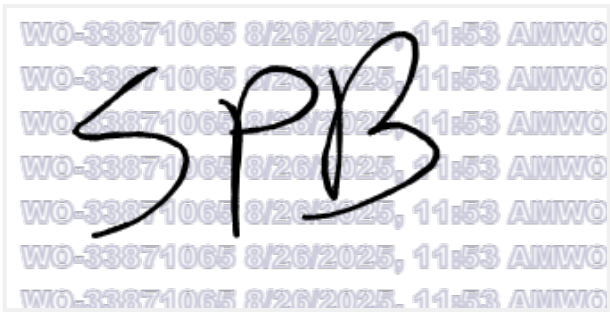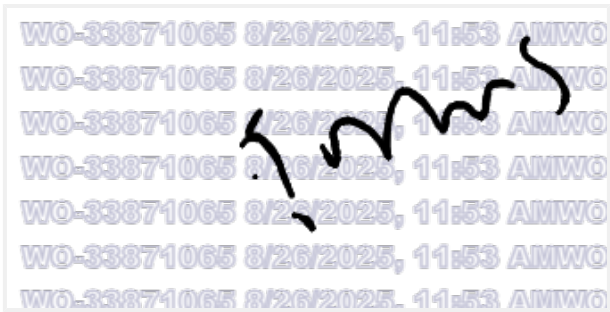

By signing this document

- Customer representative confirms they have reviewed the work carried out as described and accepts the system for use.
- Roche Support person confirms the above mentioned work and checks have been completed..
- .